



## Potato Dextrose Agar

### Intended use

Potato Dextrose Agar is recommended for the isolation and enumeration of yeasts and moulds from water, dairy, other food products and clinical samples.

### Composition\*\*

Ingredients	Gms / Litre
Potatoes, infusion from	200.000
Dextrose	20.000
Agar	15.000
Final pH ( at 25°C)	5.6±0.2

\*\*Formula adjusted, standardized to suit performance parameters

### Directions

Suspend 39.00 grams in 1000 ml purified/ distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Mix well before dispensing. In specific work, when pH 3.5 is required, acidify the medium with sterile 10% tartaric acid. The amount of acid required for 100 ml of sterile, cooled medium is approximately 1 ml. DO NOT HEAT the medium after addition of the acid.

### Principle And Interpretation

Potato Dextrose Agar is recommended by APHA (9) and F.D.A. (4) for plate counts of yeasts and moulds in the examination of foods and dairy products (3). Potato Dextrose Agar is also used for stimulating sporulation, for maintaining stock cultures of certain dermatophytes and for differentiation of typical varieties of dermatophytes on the basis of pigment production (8). It is also recommended by USP (10), BP (2) ,EP (3) and JP (6) for growth of fungi.

Potato infusion and dextrose promote luxuriant fungal growth. Adjusting the pH of the medium by tartaric acid to 3.5, inhibits the bacterial growth. Heating the medium after acidification should be avoided as it may hydrolyse the agar which can render the agar unable to solidify.

### Type of specimen

Food and dairy samples; Water samples, Clinical samples- skin scrapings

### Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (5,7).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (4,9,11).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards.(1)

After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions

In Vitro diagnostic Use. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets

### Limitations :

1. Heating the medium after acidification should be avoided as it may hydrolyse the agar which can render the agar unable to solidify.

### Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

### Quality Control

#### Appearance

Cream to yellow homogeneous free flowing powder

**Gelling**

Firm, comparable with 1.5% Agar gel

**Colour and Clarity of prepared medium**

Light amber coloured clear to slightly opalescent gel forms in Petri plates.

**Reaction**

Reaction of 3.9% w/v aqueous solution at 25°C. pH : 5.6±0.2

**pH**

5.40-5.80

**Cultural Response**

Cultural characteristics observed after an incubation at 22 - 25°C for 4 - 5 days .

Organism	Growth	Ascospore formation
<i>Aspergillus brasiliensis</i> ATCC 16404 (00053*)	luxuriant	Negative
<i>Candida albicans</i> ATCC 10231 (00054*)	luxuriant	Negative
<i>Saccharomyces cerevisiae</i> ATCC 9763 (00058*)	luxuriant	Positive

Key : \* - Corresponding WDCM numbers

**Storage and Shelf Life**

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Use before expiry date on the label.

Product performance is best if used within stated expiry period.

**Disposal**

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (5,7).

**Reference**


1. Baird R.B., Eaton A.D., and Rice E.W., (Eds.), 2015, Standard Methods for the Examination of Water and Wastewater, 23rd ed., APHA, Washington, D.C.
2. British Pharmacopoeia, 2017, The Stationery office British Pharmacopoeia
3. European Pharmacopoeia, 2017, European Dept. for the quality of Medicines.
4. FDA Bacteriological Analytical Manual, 2005, 18th Ed., AOAC, Washington, DC.
5. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2<sup>nd</sup> Edition.
6. Japanese Pharmacopoeia, 2016.
7. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
8. MacFaddin J. F., 1985, Media for the Isolation-Cultivation-Identification-Maintenance of Medical Bacteria, Vol.1, Williams and Wilkins, Baltimore
9. Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
10. The United States Pharmacopoeia, 2018, The United States Pharmacopoeial Convention. Rockville, MD.
11. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.



In vitro diagnostic medical device



CE Marking

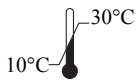


Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,  
23 Vadhani Industrial Estate,  
LBS Marg, Mumbai-86, MS, India

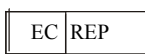

CE Partner 4U, Esdoornlaan 13, 3951  
DB Maarn The Netherlands,  
[www.cepartner4u.eu](http://www.cepartner4u.eu)

#### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.